


DAVOL* PERMASORB* DISPOSABLE FIXATION DEVICE

Absorbable Prosthetic Fixation

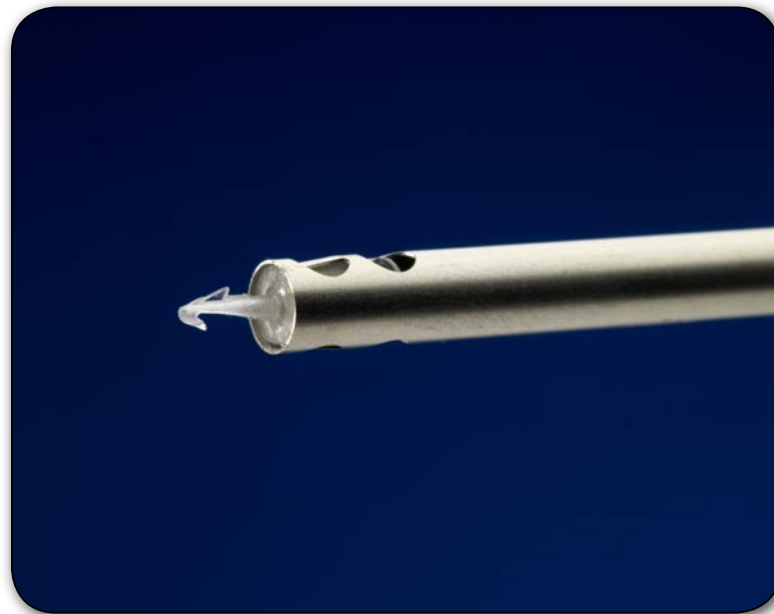

Davol* PermaSorb* Disposable Fixation Device is just one in a complete family of hernia repair products:

Ventral Hernia Repair Products

- Bard* Composix* E/X Mesh
- Bard* Composix* Kugel* Hernia Patch
- Bard* Ventralex* Patch
- Bard* Dulex* Mesh
- Bard* Composix* L/P Mesh
- Bard* **CollaMend*** Implant
- Bard* **AlloMax*** Surgical Graft

Inguinal Hernia Repair Products

- Bard* Mesh Flat Sheets and Pre-Shapes
- Bard* PerFix* Plug
- Bard* Modified Kugel* Patch
- Bard* 3DMax* Mesh
- Bard* **Soft** Mesh

Hiatal Hernia Repair Products

- Bard* CruraSoft* Patch

Fixation Products

- Davol* Salute* II Disposable Fixation Device

Specialty Products

- Bard* **C*** Parastomal Patch

Davol Surgical Services Clinical Training Program

National training centers offer instruction in surgical techniques and the ability to view live surgery.

Speaker Program

Educational presentations are given at Grand Rounds, Society Meetings and other venues.

Procedure Introduction Kits

Video programs that describe specific hernia repair techniques and their benefits to you, your patients and your surgical practice.

These services are available for many of the Davol Hernia Repair Products. Please ask your representative.



A **BAIRD** Company

Davol Inc. • Subsidiary of C. R. Bard, Inc.
P.O. Box 8500 • Cranston, RI 02920
1.800.556.6275 • www.davol.com

Benefits of the Davol* PermaSorb* Fixation Device:

- Poly (D, L) – lactide (PDLLA) fasteners fully absorb at approximately 16 months post implantation, leaving no foreign material in the patient's body
- Absorption occurs at a near constant rate for approximately 9 months
- Introducer needle effortlessly pilots fastener through mesh and tissue
- Pre-loaded, single-use device
- Procedure-specific fixation with 5-shot and 12-shot devices available in both laparoscopic and open configurations.
- Fully disposable design
- Quick and easy to use
- Safe and effective


PermaSorb*



The Davol* PermaSorb* Disposable Fixation Device is designed for prosthetic mesh fixation and the approximation of soft tissue in laparoscopic and open surgeries. The product is a pre-loaded fixation device with absorbable PDLLA fasteners that are deployed via a 5mm, low-profile delivery system.



Davol* PermaSorb* Fastener Deployment Technology

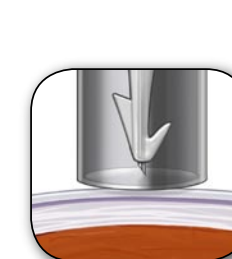

Device is supplied with pre-loaded absorbable fasteners.



Introducer needle pilots fastener through mesh and tissue.



Fastener releases from tip and securely fixates mesh into tissue.

Laparoscopic Surgical Procedures

0113090 PermaSorb* Laparoscopic Fixation Device – 5 Absorbable Fasteners

0113092 PermaSorb* Laparoscopic Fixation Device – 12 Absorbable Fasteners

Open Surgical Procedures

0113097 PermaSorb* Open Fixation Device – 5 Absorbable Fasteners

0113099 PermaSorb* Open Fixation Device – 12 Absorbable Fasteners

Please add Davol* PermaSorb* Disposable Fixation Device to my preference card.

I would like to have Davol* PermaSorb* Disposable Fixation Device in stock.

Surgeon's Signature _____

Purchase Order Number _____

Catalog Number _____

Date _____ Quantity _____


* Bard, Composix, Kugel and PermaSorb are trademarks and/or registered trademarks of C.R. Bard Inc. or an affiliate. Consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and instructions for use.

Copyright © 2007 Davol Inc. All rights reserved.
MMPSSS




A **BAIRD** Company

Davol* PermaSorb* Disposable Fixation Device with Bard* Composix* L/P Mesh in Laparoscopic Ventral Procedures

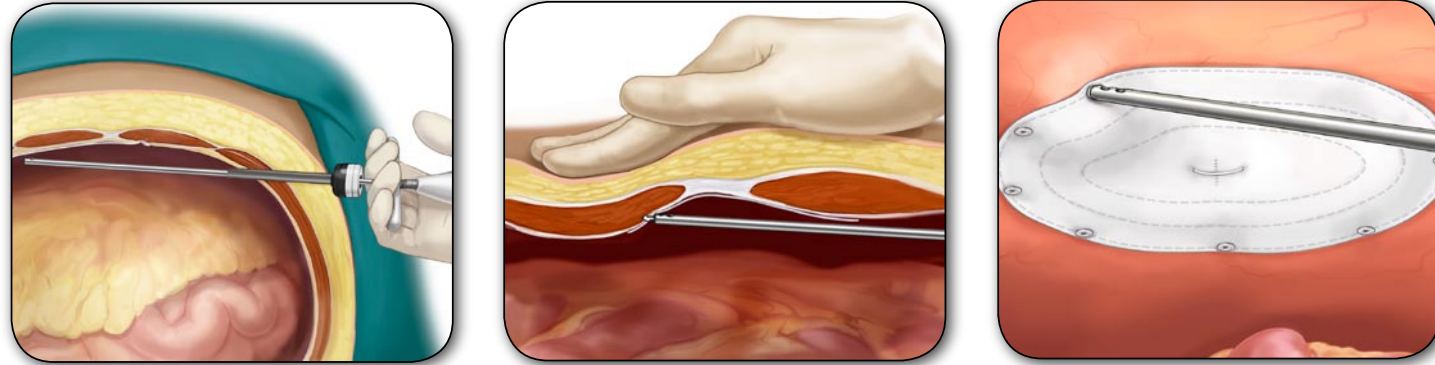

- The PermaSorb* Disposable Fixation Device provides secure fixation in laparoscopic ventral procedures
- PDLLA fasteners are fully absorbable
- Can be used with most 5 mm trocars


Davol* PermaSorb* Disposable Fixation Device with Bard* Soft Mesh in Laparoscopic Inguinal Procedures

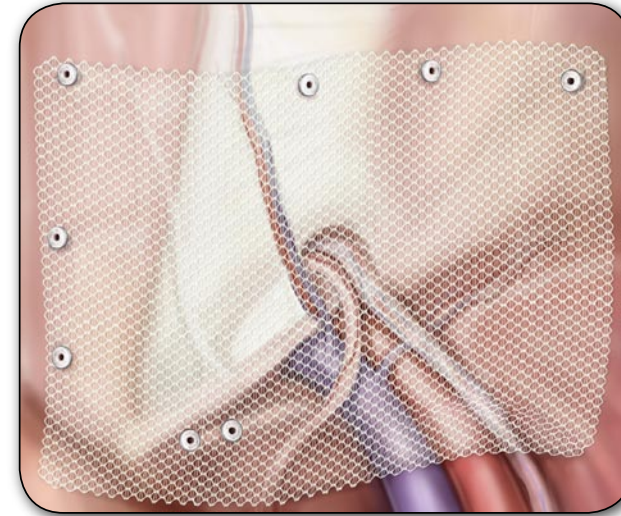


- The PermaSorb* Disposable Fixation Device provides controlled and secure fixation in laparoscopic inguinal procedures
- Provides the appropriate count of fasteners for the procedure being performed

Davol* PermaSorb* Disposable Fixation Device with the Bard* Composix* Kugel* Hernia Patch in Open Ventral Procedures

- Use of the PermaSorb* Disposable Fixation Device instead of suturing on the periphery increases procedure speeds
- Extensive lateral dissection is not required when using the Bard* Composix* Kugel* patch and the PermaSorb* Disposable Fixation Device

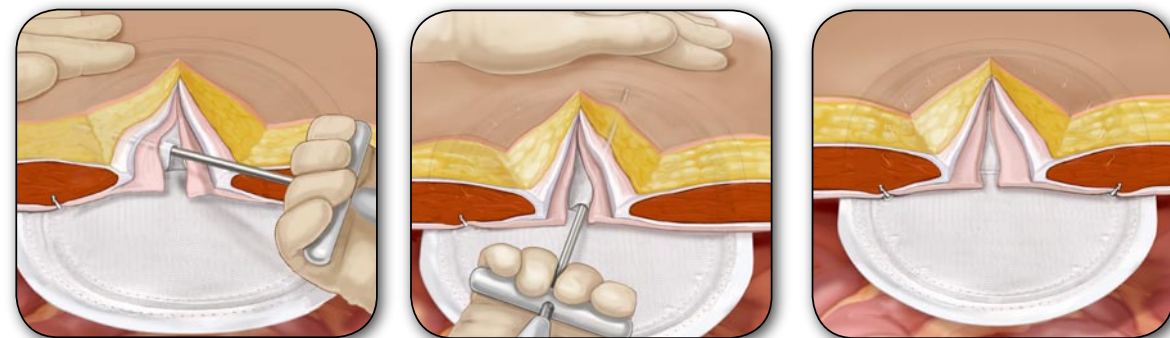

Davol* PermaSorb* Disposable Fixation Device with Bard* Soft Mesh in Open Inguinal Procedures

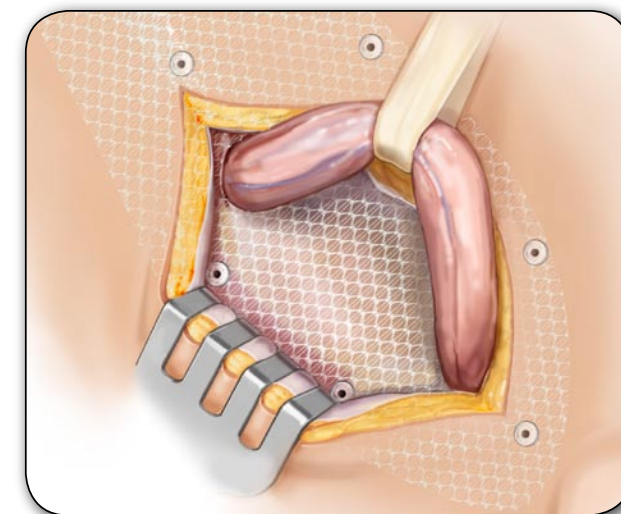

- Use of the PermaSorb* absorbable fasteners instead of suturing reduces OR time
- PDLLA fasteners are absorbable over time, leaving no foreign material in the body